


fitAims

Skincare platform Hydrificent

4,600 EUR

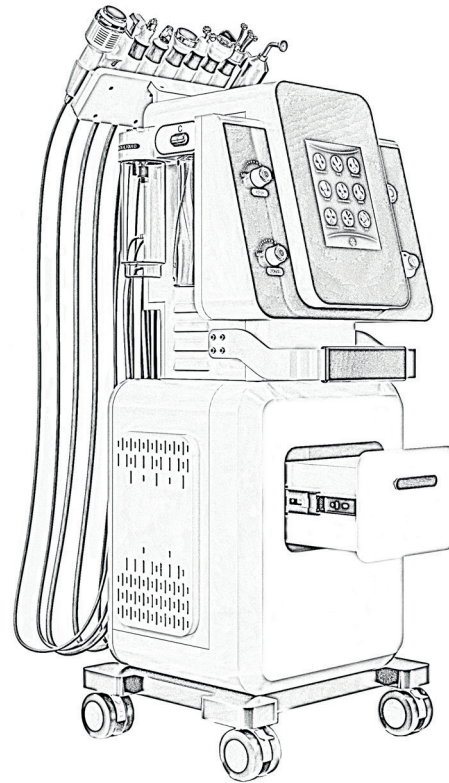

**Multifunctional skincare
platform for deep
cleaning, skin tightening
and skin rejuvenation**

8-in-1 skincare solution for the most
common skin issues

CE
certified

+421 353 211 408

info@fitaims.com

www.fitaims.com



An all-in-one professional beauty platform for your beauty studio is an optimal choice when you are torn between the desire to expand the list of services and uncertainty as to the customer demand for this or that treatment.

With 8 most popular technologies at hand you can be sure that at least 5 or 6 of them will be able to create a steady client flow. Moreover, multifunctional solutions allow saving costs and space.

What are the benefits of Hydraficient platform?

Excellently chosen treatments which can be smoothly switched from one to another. On the one hand, an operator is free to create the best treatment protocol, on the other, a customer is more likely to add an extra treatment if it can be performed almost immediately after the end of his planned procedure.

Cost-efficiency. The cost of each individual beauty device is by 20-30% higher than the same function as a part of a cosmetological combine. Therefore the investment return period will also be shorter.

Profitability. Non-invasive nature, low cost and superb results make beauty treatments top-rated among the customers. These procedures with a relatively low entry threshold give you an opportunity to offer your customers various other treatments in future.

What treatments will you be able to offer?

Hydrodermabrasion.

Safe and effective treatment for gentle exfoliation and hydration of the facial skin. An alternative of microdermabrasion with immediate results and no side effects. A hydro dermabrasion handle is equipped with a variety of treatment heads to work with different skin conditions.

Miscrocurrent facial toning.

Uses weak electrical signals to stimulate the muscles and skin. Promotes wound healing, reduces inflammation, boosts circulation, and improves muscle function.

Ultrasonic skin scrubbing.

Handle-emitted soundwaves help loosen and unclog the pores. The ultrasonic scrubbing spatula slides over the skin, removing blackheads, softening skin keratinocytes, and deeply cleansing skin dirt.

Ultrasonic phonophoresis.

An effective and painless treatment, making cosmetic products easier to penetrate into the stratum corneum. Used for intense skin hydration, saturation of the skin with minerals, vitamins and hyaluronic acid. It makes the skin smooth and supple, returns a healthy look, tone and removes fine wrinkles.

Deep skin hydration with atomized solution.

A beauty preparation is converted into a mist-like solution through high-pressure atomization and then quickly poured into the dermis by high-pressure spraying to achieve a deep hydrating effect.

Deep skin cleaning with a bubble mousse handle.

The treatment aimed at skin cleaning, moisturizing and rejuvenation with the aid of a dense foam with a large amount of active oxygen and hydroxide ions penetrating into the skin.

RF skin tightening.

RF treatment is a non-invasive effective anti-aging procedure. Involves deep skin layers heating to promote collagen production, to smooth and brighten the skin. The results of a course of 6-10 treatments may last up to 12 months.

Cryo-thermo-lifting.

Thermal massage with a special cold/hot hammer. Cold Hammer 0°C: skin calming and toning, shrinking pores. Warm Hammer 45°C: heating the skin surface, accelerating blood circulation, and increasing the absorption rate when used with beauty products.

Highlights



8-in-1 technology
Well-chosen highly-combinable functions



Intelligent user interface
Preset treatment parameters for flawless operation



8.4" screen
Comfortable operating experience



Compact design
Space-saving solution for your beauty studio



Scope of treatment



Deep skin cleaning

Removing blackheads, grease, dust, dead skin cells, cleaning pores, removing mites



Facial skin rejuvenation

Skin moisturizing, shrinking pores, skin whitening, acne treatment



Skin tightening

Face lifting, eye lifting, wrinkle removal, overall improvement in skin structure and texture due to boosted collagen

Accessories



Cold & Warm Hammer

120 EUR



Ultrasound scrubbing handle

120 EUR



Ultrasonic phonophoresis handle

120 EUR



Atomizing handle

120 EUR



Multipole microcurrent handle

120 EUR



RF skin tightening handle

120 EUR



Two-pole microcurrent handle

120 EUR



Hydrodermabrasion handle

120 EUR



Bubble mousse handle

120 EUR



Dermabrasion heads (10 pcs) for hydrodermabrasion handle

30 EUR

Consumables

Cosmetic solutions
Purified water for diluting the solutions

Downtime

12 hours of continuous work guaranteed

Warranty

2-year warranty
During the first year of warranty, spare parts and delivery are free of charge.
During the second year, only spares are free of charge, the delivery should be covered by the customer.

Servicing/Maintenance

The machine is maintenance-free.
Keep the handles and bottles clean to ensure hygiene.

Specifications

RF frequency	3Mhz
Vacuum	0 - 90 kPa
Flow capacity	30 L/max
Power requirements	100V-240V 50Hz-60Hz
Number of dermabrasion heads	10
Screen size	8.4"
Power consumption	600 W
Dimensions	60*50*180 cm
Packing details	54*57*115 cm
Net weight	65 kg

